
# Vertidisan®

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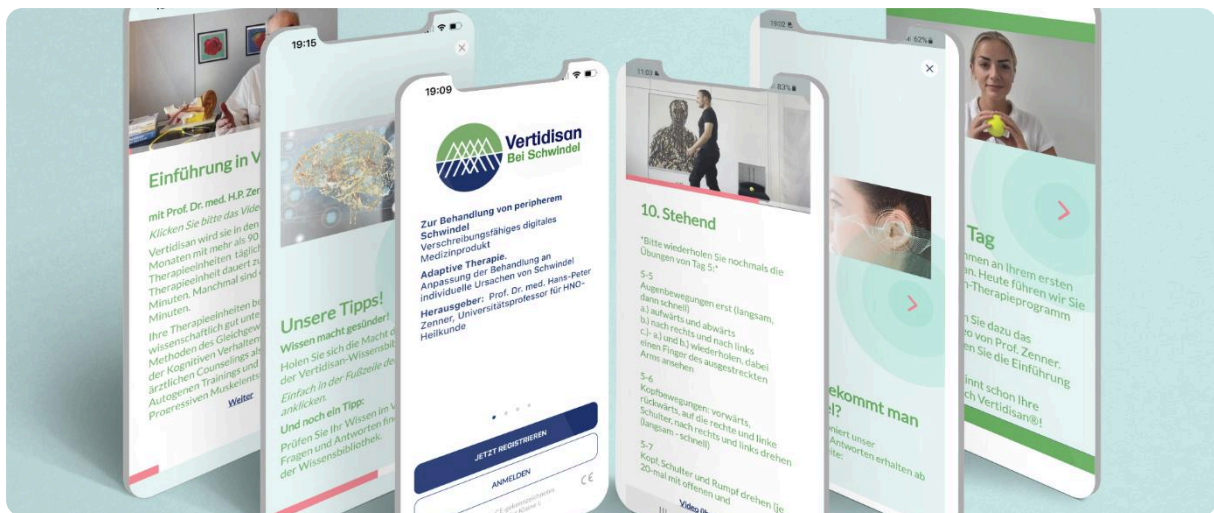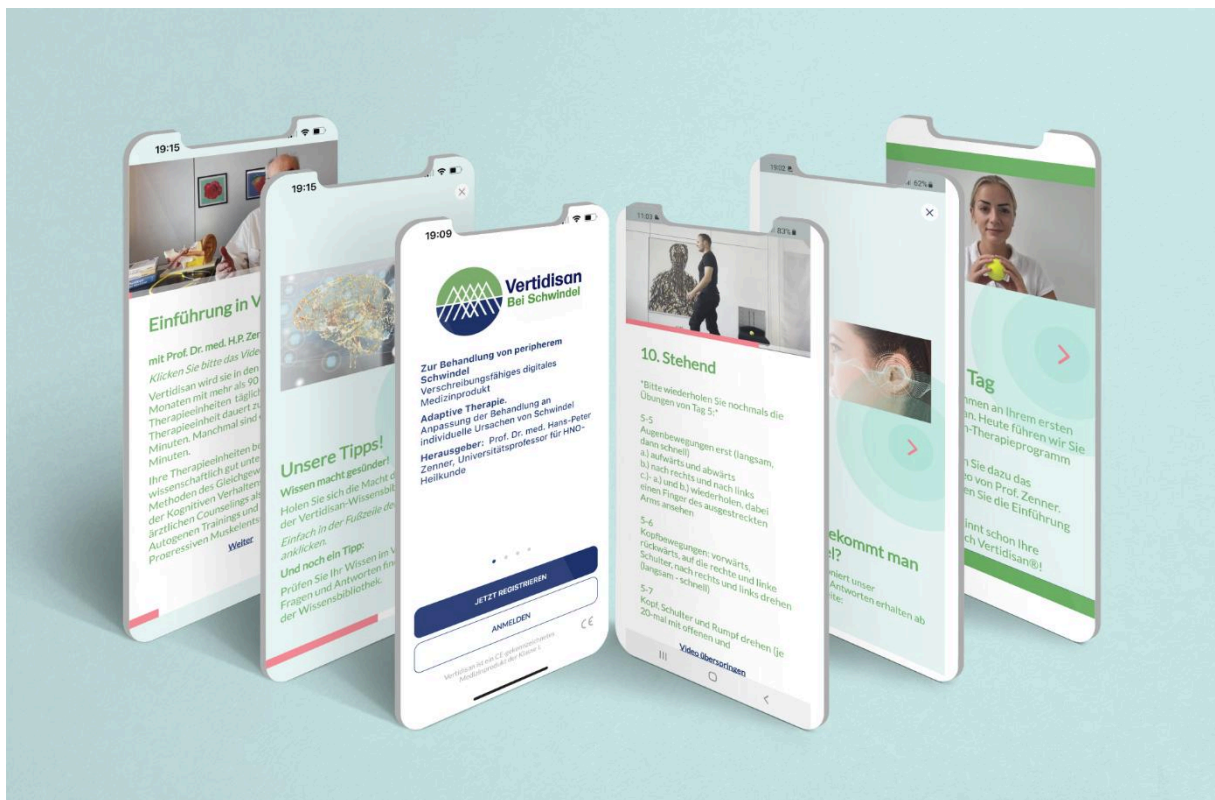
Digitineers GmbH & Co KG Vertidisan Vertigo App Instructions for Use

## Table of Contents

|                                                                      |    |
|----------------------------------------------------------------------|----|
| <b>Introduction</b>                                                  | 5  |
| <b>1. Notes on Use</b>                                               | 5  |
| Intended Purpose                                                     | 5  |
| Medical Purpose                                                      | 5  |
| Safety Information                                                   | 6  |
| Additional Information                                               | 6  |
| <b>2. Technical Requirements</b>                                     | 7  |
| <b>3. How It Works</b>                                               | 7  |
| <b>3.1 The Contents</b>                                              | 7  |
| <b>3.2 Functions</b>                                                 | 8  |
| <b>3.3 Product Description</b>                                       | 8  |
| Functionalities                                                      | 9  |
| <b>3.4. Introduction and Programme Configuration</b>                 | 9  |
| <b>3.5. Daily Therapy</b>                                            | 10 |
| <b>3.6. Progress Log and Sharing of Progress Data</b>                | 10 |
| <b>3.7. Scientific Data Evaluation</b>                               | 10 |
| <b>4. Legal Requirements for the Use and Processing of Your Data</b> | 10 |
| <b>5. Medical Requirements</b>                                       | 11 |
| <b>5.1. Medical Requirements for Use</b>                             | 11 |
| <b>5.2 Contraindications, Warnings, Interactions</b>                 | 11 |
| Contraindications                                                    | 11 |
| Pregnancy                                                            | 12 |
| Driving and Operating Machinery                                      | 12 |
| Warnings                                                             | 12 |
| Interactions                                                         | 12 |
| <b>5.3. Accessibility</b>                                            | 12 |
| <b>6. Risks and Side Effects</b>                                     | 13 |
| <b>6.1. Risks</b>                                                    | 13 |
| <b>6.2. Side Effects</b>                                             | 14 |

|                                                                                                  |           |
|--------------------------------------------------------------------------------------------------|-----------|
| Known Hazards and Side Effects .....                                                             | 14        |
| Other Hazards or Side Effects .....                                                              | 15        |
| Reporting a Risk or Side Effect .....                                                            | 15        |
| <b>7.0 Installing Your Vertidisan App .....</b>                                                  | <b>16</b> |
| <b>8.0. Operating Vertidisan .....</b>                                                           | <b>17</b> |
| <b>8.1 First Steps from Opening the App to the Start of Therapy (so-called Onboarding) .....</b> | <b>17</b> |
| <b>8.2 Start of Therapy .....</b>                                                                | <b>18</b> |
| <b>8.3 Continuation of Therapy .....</b>                                                         | <b>19</b> |
| <b>8.4 End of Therapy .....</b>                                                                  | <b>19</b> |
| <b>8.5 Controlling the App .....</b>                                                             | <b>19</b> |
| Next Button, Forward and Back Arrows .....                                                       | 19        |
| Footer Bar/Footer .....                                                                          | 19        |
| Knowledge Icon in the Footer .....                                                               | 20        |
| Menu Icon in the Footer .....                                                                    | 21        |
| <b>8.6 Intended Use .....</b>                                                                    | <b>22</b> |
| <b>9.0. Classification and Operating System Requirements .....</b>                               | <b>23</b> |
| <b>10.0. Performance Characteristics .....</b>                                                   | <b>23</b> |
| <b>11.0 Compatibility Commitments Regarding Supported Platforms .....</b>                        | <b>23</b> |
| <b>12.0 CE Marking .....</b>                                                                     | <b>23</b> |
| <b>13.0 Data Protection .....</b>                                                                | <b>23</b> |
| <b>14.0 Product Labelling .....</b>                                                              | <b>24</b> |
| <b>15.0 Contact .....</b>                                                                        | <b>24</b> |
| <b>16.0 Manufacturer .....</b>                                                                   | <b>24</b> |
| <b>Annex 1: .....</b>                                                                            | <b>25</b> |
| <b>Indications .....</b>                                                                         | <b>25</b> |
| Vertigo Due to Disorders of Vestibular Function .....                                            | 25        |
| <b>Contraindications .....</b>                                                                   | <b>25</b> |

## Introduction

These instructions for use are intended to ensure the correct use of Vertidisan®. Please read this information carefully and follow the instructions it contains. In addition, the Vertidisan® support team is available for enquiries at support@digitineers.de during normal business hours.

## 1. Notes on Use

### Intended Purpose

Vertidisan® serves to improve vertigo in peripheral vestibular disorders. Vertidisan also serves to document manually entered health data and to forward this documentation.

Vertidisan is not intended to produce or provide information that is used to make decisions for diagnostic or therapeutic purposes.

### Medical Purpose

Vertidisan® ("App") is an online-based, CE-marked medical device that supports the therapy of patients with vertigo resulting from a peripheral vestibular disorder. Vertidisan® contains digital, interactive, multimodal and diagnosis-adapted therapeutic training programmes consisting of 90 therapy sessions with vertigo exercises, behavioural therapy units, vertigo counselling and instructions for muscle relaxation and muscle strengthening. These are complemented by relaxation phases and health education. Vertidisan® is suitable for persons affected by so-called peripheral vertigo, vestibular neuronitis, otogenic vertigo, aural vertigo, other disorders of vestibular function and vertiginous syndromes NOS (not otherwise specified). A detailed catalogue of indications and contraindications for the App can be found in **Annex 1** of these instructions.

The training can be carried out without supervision by a doctor or therapist. Before use, however, a doctor should have confirmed that Vertidisan® is suitable for you and your health situation. Unless otherwise instructed by a doctor, one Vertidisan® training unit should be carried out daily. The duration of use depends on the type of complaint and/or the duration of the illness and is therefore determined solely by your doctor. The App supports the treatment and relief of vertigo both within the framework of care provided by a doctor and when used independently without a doctor. However, it does not provide any medical or therapeutic diagnosis, recommendation or medical treatment.

The therapy programme of this App is edited by Prof. Dr. med. Hans-Peter Zenner, University Professor of Otorhinolaryngology (ENT), Tübingen. The effectiveness of the App has been scientifically tested; nevertheless, it does not offer the best therapy in every case. Vertigo can also be a sign of serious illness. Some patients are unable to follow the exercise programme.

### **Minimum Age**

18 years.

**Please read these instructions carefully before using the App. In addition, a doctor should have confirmed to you before use that Vertidisan® is suitable for you and your particular health situation.**

Questions going beyond these instructions can be sent at any time to [support@digitineers.de](mailto:support@digitineers.de).

### **Safety Information**

- In individual cases, vertigo may be due to serious illness. We therefore recommend seeking medical advice before using Vertidisan.
- Individual exercises are not suitable for every patient. We therefore recommend seeking medical advice before using Vertidisan.
- Clarify with your doctor whether performing physical exercises on your own is suitable for you.
- Do not exercise under the influence of alcohol or narcotics.
- Choose a suitable location with sufficient space.
  - Avoid tripping hazards and unsuitable clothing.
  - Have something to hold on to (e.g. a chair) within easy reach.

### **Additional Information**

- During the exercises, please put your smartphone down (e.g. on a stand).
- During all exercises, the App guides you acoustically so that you do not have to look at the screen constantly.
- During and after the first exercises, temporary vertigo and/or muscle soreness may occur – this is normal at the beginning and disappears later.
- The exercises usually take approx. 10–15 minutes a day, sometimes a little longer.

- Medical behavioural therapy and medical counselling can take up to 30 minutes. You may therefore wish to schedule your training days at the weekend, for example.
- In the App you can watch the instruction videos for earlier exercises as often as you like by reopening the tiles of days whose programme you have already completed.
- In most cases, you should carry out the exercises for 3 months. If you are free of vertigo earlier, you can click the “End of therapy” button (in the library) at any time and then stop.
- From the age of 60 in particular, more than 3 months of therapy are required.
- If you have consented to the use of your health data, the App synchronises the usage data with the server. You can withdraw your consent to the use of your health data at any time. The App will continue to function. However, all your personal health data will be irrevocably deleted.
- If you have consented to the use of your health data, your data can be viewed again if you change your smartphone/tablet.
- The App is not designed for simultaneous use by several users or on several devices.

## **2. Technical Requirements**

Vertidisan® can be accessed as an interactive online programme via any internet connection. To use Vertidisan®, the patient receives an access key. This access key entitles the patient to use Vertidisan® for 90 days of use after initial registration.

You need a smartphone or a tablet PC (smartphone or tablet with Android operating system version 12 or higher; iOS operating system version 18 or higher). You also need an e-mail address. As Vertidisan® guides you through the programme, no special computer skills are required.

## **3. How It Works**

### **3.1 The Contents**

The medical contents of Vertidisan® were compiled by experienced doctors. All contents and treatment principles on which Vertidisan® therapy is based are derived from science, from treatment guidelines for psychotherapists and doctors or from recommendations of medical societies for the treatment of vertigo. In the treatment units, most of the therapy units are prepared by a doctor or therapist as audio or video recordings.

## 3.2 Functions

Vertidisan® comprises recognised treatment elements of physiotherapy, cognitive behavioural therapy, counselling and muscle relaxation as well as muscle strengthening. The App provides information on diseases involving vertigo, presents exercises, methods and techniques and supports you with additional written information to help you consolidate the contents of the Vertidisan® exercises and integrate them into everyday life.

Vertidisan® has several main and secondary functions. Vertidisan® offers the user numerous videos, each with supplementary text information, providing instructions for therapeutically effective exercises for vertigo. Vertidisan® then provides the next instruction, which the user can in turn learn and carry out, and so on. This creates a communication flow in which the other functions are additionally embedded. The second main function consists of therapeutically supportive video sequences with units of medical cognitive behavioural therapy. The third and fourth main functions are medical vertigo counselling and muscle relaxation by means of Jacobson's progressive muscle relaxation, autogenic training or short relaxation.

Further functions are available when using Vertidisan®. These include:

- Audios, e.g. with relaxation instructions
- A health knowledge library
- Recording of the therapy outcome status (voluntary)
- Summaries of therapy progress as a progress report for self-monitoring (as a PDF) and/or in machine-readable form for digital transfer to the doctor's or therapist's digital practice system as well as to the electronic patient record (voluntary export function).

The user can interrupt Vertidisan® at any time and then continue at the point where he/she left the programme. Vertidisan® is password-protected and available via the internet at any time. This allows use independent of time, place and device. The statutory data protection requirements are met. As a German manufacturer, we are subject to the German Federal Data Protection Act (Bundesdatenschutzgesetz). The applicable requirements under this law are met.

## 3.3 Product Description

Vertidisan® offers a structured therapy for patients with peripheral vestibular disorders. The App has no measuring or diagnostic functions and is therefore a Class I medical device under the Medical Device Regulation (European Medical Device Regulation).

In accordance with German scientific guidelines and international publications, the therapy programme contains therapy units from five main areas:

- **Vertigo exercises:** At the heart of the App are video demonstrations of scientifically based vertigo exercises as an almost daily item of the therapy programme. The patient follows the video demonstration by a therapist and additionally receives written instructions.
- **Cognitive behavioural therapy:** An effective component of the App are units of cognitive behavioural therapy, which are completed as a manualised, structured CBT programme.
- **Vertigo counselling:** In medical video sequences and short articles, users receive information about vertigo and about current treatment and diagnostic options.
- **Muscle relaxation:** This part contains an audio file with medical instructions for Jacobson's progressive muscle relaxation and autogenic training.
- **Otago Exercise Program (OEP):** In patients with a tendency to fall, the "Otago Exercise Program" carefully strengthens the muscles, particularly of the legs. By treating not only balance but also the musculature in this way, the tendency to fall is further reduced.

Offered in addition to the therapy:

- **Knowledge library** – A knowledge library on many questions relating to vertigo is also available.

## Functionalities

- Introduction to configure the App. Confirmation by the user that, from a medical point of view, he/she is able to exercise independently.
- Management of consents
- Daily therapy programme with therapy units
- Video or audio instructions for the therapy units mentioned above
- Text instructions for the therapy units mentioned above
- Access to the health knowledge library
- Input for the therapy log and output of the log
- Export and share. Data on the progress achieved can be exported and shared with a doctor/therapist at any time.
- Access to information required by regulation. This is accessible via the menu.

## 3.4. Introduction and Programme Configuration

In the introduction to the App, you can answer questions about yourself, your complaints and any existing medical findings, among other things. There is no

obligation to answer the health-related questions. You can also use the App without providing this information.

### 3.5. Daily Therapy

After setting up your Vertidisan® therapy account, you will regularly receive instructions for your daily therapy in the App.

**Stop an exercise if it triggers or worsens vertigo or pain, or if you notice other side effects. See a doctor if these symptoms persist.**

### 3.6. Progress Log and Sharing of Progress Data

If you wish, you can keep a therapy outcome log with the App. To do so, you can enter information about your vertigo. You can display the entries in your therapy log as a progress report in PDF format at any time. You can then make the PDF file available to your doctor or therapist.

### 3.7. Scientific Data Evaluation

You can provide us with the data of the therapy report in anonymised form, i.e. without any reference to your person, so that a scientific evaluation to improve the App is possible. Anonymous means that the data contain neither your name, nor your date of birth, nor your address, nor any code linked to you.

## 4. Legal Requirements for the Use and Processing of Your Data

As a prerequisite for use, you confirm at the start of the introduction that you

- are at least 18 years old
- have read and accept the Vertidisan® General Terms and Conditions (terms of use)
- have read and accepted the Vertidisan® privacy policy.

It is not a prerequisite that you consent to the processing of your personal health data by Vertidisan®. As the manufacturer, we generally have no access to your personal health data.

If you wish the App to record and evaluate your health-related data for you, you must first give your consent. The purpose of the data processing is to help you document your treatment data for yourself and, if you wish, also for your doctor.

You can withdraw your consent to the processing of your data at any time and end the processing.

## 5. Medical Requirements

### 5.1. Medical Requirements for Use

Vertidisan® (“App”) is an online-based, CE-marked medical device that supports the therapy of patients with vertigo resulting from a peripheral vestibular disorder.

**If you experience pain, increased vertigo, malaise or vomiting during or after using the App, you should stop the programme and, if necessary – at your own discretion – consult a doctor.**

Pause Vertidisan® in the case of:

- fever
- general feeling of illness

### 5.2 Contraindications, Warnings, Interactions

This section contains several very important sub-items:

#### Contraindications

Contraindications are the diseases or complaints in which Vertidisan® must not be used.

- Diseases of the cervical spine and neck
- Acute inflammatory diseases
- Febrile condition
- Heart failure (cardiovascular decompensation)
- Lung failure (pulmonary insufficiency)
- Motor disorders
- Seizure disorders (e.g. epilepsy)
- Spinal cord damage
- Psychoses
- Skin sensitivity disorders with walking impairment
- Pregnancy
- Thrombosis

- Blood coagulation disorders, including anticoagulant therapy
- Large skin defects
- Burns requiring hospitalisation
- Metastases of malignant tumours
- Severe visual impairment
- Frailty, frailty syndrome

## Pregnancy

Vertidisan® must not be used in the event of pregnancy.

## Driving and Operating Machinery

As long as Vertidisan® does not trigger vertigo, the App does not impair your ability to drive or operate machinery. Vertigo can impair your ability to drive or operate machinery.

## Warnings

Observe the warnings. They refer to complaints or illnesses that it is important for you to have discussed with your treating doctor, because special attention is required in these cases. You will find the warnings in the instructions for use.

## Interactions

No interactions of Vertidisan® with additionally used medicines, foods or food supplements are known. Nevertheless, Vertidisan® should not be used together with medicines whose intake is associated with a risk of falling. Ask your doctor or pharmacist about this.

## 5.3. Accessibility

Up-to-date information on accessibility can be found in the knowledge library in the section “Accessibility”.

**The App is not fully accessible; specifically, it is not suitable for users with severe visual impairments or with severe motor disorders.**

For people with **hearing impairments**, all spoken texts are also available in written form. If you have a suitable hearing aid, middle ear implant or cochlear implant, you can stream the spoken texts directly to your device or implant.

For people with **visual impairments**, we have designed the control buttons to be large. You can also listen to the text of all therapy units.

For people with **motor disorders of the hand**, we have designed the control buttons to be sufficiently large.

## 6. Risks and Side Effects

It is possible that not every patient will benefit from using the programme. If working with the programme is too stressful for you, you should take a break, suspend use of the programme and talk to your doctor about it. Not every exercise is equally suitable for every patient. If an exercise causes you difficulties, feels unpleasant or impairs you, do not continue with this exercise. Consult your doctor if necessary.

### 6.1. Risks

Hazards can arise from incorrect performance of exercises.

**Caution!** Therefore, make sure that you perform the exercises as shown to you in each case. If this is not possible for you, do not perform the respective exercise.

**Caution!** Existing medication and/or treatments must not be changed without consulting your doctor/psychotherapist.

**Caution!** Always follow your doctor's instructions regarding the exercises you are allowed to perform in your particular state of health. Your doctor's instructions always take precedence over the instructions in the App and must always be followed exactly. Disregarding your doctor's instructions when using Vertidisan® can lead to pain, vertigo and injuries.

**Caution!** Before performing an exercise, always pay close attention to the instructions in the instruction video and always read the explanatory texts accompanying the instruction videos. They contain important information on the correct performance of the exercises. If you are not sure how to perform an exercise, stop and ask your doctor. Although the exercises have been selected so that they are suitable for independent practice at home without a doctor present, incorrect performance of the exercises can lead to pain, vertigo and injuries.

## **6.2. Side Effects**

All therapeutic procedures can have side effects; however, these do not occur in every patient.

Before this App was used for the treatment of diseases and complaints, clinical evaluations and applications regarding the effectiveness, hazards and side effects of the App were carried out. The latter may be “known”, or they may not yet have occurred but be foreseeable for doctors. These are referred to as “other hazards and side effects”. They are listed after the known hazards and side effects:

### **Known Hazards and Side Effects**

- Temporary vertigo, increase in vertigo or balance disorder. An exercise may initially trigger vertigo or unsteadiness of gait, which, however, soon subsides
- Temporary nausea, loss of appetite
- Temporary tiredness
- Temporary headache
- Temporary neck pain
- Temporary blurred or tunnel vision
- Temporary pain in the musculoskeletal system
- Tension in the muscles or back
- Temporary tiredness
- Temporary knee pain
- Temporary tinnitus
- Occasionally longer-lasting, temporary increase in vertigo

**Caution!** If the App triggers rotational vertigo or unsteadiness of gait that does not subside, please consult your doctor.

### **Other Hazards or Side Effects**

The following hazards and side effects are foreseeable but have not yet been observed:

#### ***Hazards***

Incorrect performance of exercises due to

- Font too small
- Confusing description
- Unclear handling/operation

#### ***Side effects***

The following side effects are possible:

- Cardiovascular syncope
- Increase in blood pressure
- Connective tissue injuries in the musculoskeletal system, muscle injuries
- Oedema (swelling)
- Inflammation of the nasal mucosa (rhinitis)
- Aversion to psychotherapy

### **Reporting a Risk or Side Effect**

You can report side effects to the manufacturer as follows:

1. By telephone on 07071-704 3011 (international: +49 7071 704 3011)
2. By e-mail to [meldung@digitineers.de](mailto:meldung@digitineers.de) (not data-secure)
3. By post (also anonymously):

Digitineers GmbH & CoKG  
August-Bebel-Str. 9  
72072 Tübingen  
Germany

Serious side effects caused by Vertidisan must be reported to the German Federal Institute for Drugs and Medical Devices (Bundesinstitut für Arzneimittel und Medizinprodukte, BfArM; [mp-vigilanz@bfarm.de](mailto:mp-vigilanz@bfarm.de)) using the reporting form <https://www2.bfarm.de/medprod/mpsv/>.

## 7.0 Installing Your Vertidisan App

You need an internet connection to install and use the App. Before installing the App, please make sure that you are either connected to a Wi-Fi network or have a large data volume with your network provider.

1. If you own an Apple device, i.e. an iPhone or iPad, you have an iOS operating system. First download the App from the App Store:

- a. To do so, open your App Store and tap the search icon.
- b. Enter "Vertidisan" in the search field.
- c. Tap the download icon.
- d. Wait until the blue circle next to Vertidisan® has finished loading, and then tap "Open".

2. If you have a smartphone or tablet with Android, download the App from the Google Play Store.

- a. To do so, open your Play Store.
- b. Enter "Vertidisan" in the search field.
- c. Tap "Vertidisan" next to the Vertidisan icon.
- d. Tap "Install".
- e. Confirm access by tapping "Accept" and open the App.

Note: All information on access permissions can be found in the Data Security section.

3. You can now open the Vertidisan App at any time by tapping the Vertidisan icon/symbol in your menu.

## 8.0. Operating Vertidisan

### 8.1 First Steps from Opening the App to the Start of Therapy (so-called Onboarding)

Open the App via the Vertidisan icon on your smartphone/tablet. Register and enter your access code.

- **Access code**

Provided the requirements are met, persons with statutory health insurance (in Germany) receive the access code free of charge. Private patients can purchase the code and submit the invoice to their insurance company for reimbursement. Details can be found at [www.verti-disan.de](http://www.verti-disan.de).

- **Consent to the terms of use**

You can read the terms of use and must give your consent.

- **Accept data protection**

You can read the data protection information and must accept it.

- **Consent to the processing of your personal data**

Here you can read information on the processing of personal data (without health data) pursuant to § 4 DiGAV (German Digital Health Applications Ordinance), and you must give your consent.

- **Voluntary consent to the processing of your health data**

Next, you will be asked for your voluntary consent to the processing of health data. With your consent, we can offer you individual therapy progress reports, which you can read yourself or give to your doctor, or transmit digitally to your doctor or to your electronic patient record (ePA). To give your consent, tap the "Confirmation" button.

*Note: You do not have to give this voluntary consent. The App also works without this consent. You can also withdraw your consent at any time in the menu under Consents. Your personal data will not be passed on to third parties, will be stored only on German servers under the highest security standards and will be used exclusively for the use and further development of the App. Further information can be found in our detailed privacy policy.*

- **Scientific questionnaires/therapy progress report**

At the push of a button, Vertidisan can create a therapy progress report for you at any time. You can read it yourself, print it out and/or forward it to your doctor as a printout or digitally. For this purpose, Vertidisan has an interface to the

usual practice management system of statutory health insurance physicians and to your electronic patient record (ePA).

The prerequisite is that you complete one or more scientific questionnaires of your choice. Select one or more of the (scientific) questionnaires and complete them (provided you have consented to the use of your health data). The time indicated for completing each questionnaire will help you choose. You can select the respective answers, where applicable, by tapping the buttons.

*Note: The more questionnaires you select and the longer the questionnaires you select, the more detailed your progress report will be. If you do not complete any questionnaire, no progress report can be created. The App will still work.*

**From day 2 onwards, tap the “Reports” button in the footer and then go to your questionnaires via the “Answer questionnaire” button.**

- **Creating your therapy progress report**

Tap “Reports” in the footer.

**Would you like to learn more about the evaluation?** Then open the tile “Therapy report: Evaluation” in the knowledge library (“Knowledge” icon in the footer) or download the PDF file there.

- **Diagnosis and selection of your therapy programme**

Finally, Vertidisan asks you for your diagnosis, which you can select (if you do not know it, select “unknown”). Depending on the diagnosis, suitable exercises for your therapeutic training as well as further therapeutic units have been compiled in the App.

- **Reminder of your exercises**

During onboarding, Vertidisan asks you whether you would like to be reminded of the daily therapy sessions by a (push) message on your smartphone/tablet. If you have given your consent, you will automatically receive a (push) message if you have not yet carried out your therapy. You can change your consent to notifications at any time in the menu.

**On your smartphone/tablet, go to “Settings” and there, for example, to “Notifications” and “Sounds/Vibration” to set the type of notification you would like.**

## **8.2 Start of Therapy**

The App now automatically takes you to the therapy start page of your therapy programme. Your therapy programme is based on the health information you

provided during onboarding. Tap the tile “Day 1”. Only on day 1 will you first receive a series of information. The therapy follows afterwards.

**For the therapy, open the respective video (or audio) and follow the video or audio.**

If you tap “Next” after the therapy, you will reach an overview of the tiles of the therapy days you have completed and of the next therapy days. You can in principle reopen the old tiles at any time; the new ones only from the following day.

In addition, a footer/footer bar becomes visible, which is described in more detail in section 8.5.

### **8.3 Continuation of Therapy**

The contents of the App build on one another. After all contents of a daily unit have been completed, the application automatically unlocks the next daily unit from midnight. To access it, click the “Today” icon in the footer (see section 8.5 Controls below).

### **8.4 End of Therapy**

You can reach the therapy unit “End of therapy” via the “Knowledge” icon (icon with a book).

**If you wish to end the therapy early, we ask you to complete the therapy unit “End of therapy” beforehand.**

### **8.5 Controlling the App**

#### **Next Button, Forward and Back Arrows**

After the therapy start page has opened and within all further therapy units, the easiest way to control the App is via the respective “Next” button or the forward and back arrows until the end of a daily programme has been reached.

**At the end of a daily programme, the footer bar described below with four icons becomes visible:**

#### **Footer Bar/Footer**

When you have completed a daily programme, a footer bar/footer appears at the bottom. In the footer bar/footer you will find the following icons:

- Today icon
- Reports icon

- Knowledge icon
- Menu icon

### **Today Icon in the Footer**

Via the “Today” icon in the footer bar you can call up the programme of the current day, of all days already completed and, from midnight, the next new daily programme via the respective day tile.

- **Tiles**

By clicking on a tile, you can call up the programme of the respective day.

- **PDF**

You can reach the PDF documents via the tiles. Some of the day tiles have a downward arrow. If you click on the arrow, you can download the PDF document.

### **Reports Icon in the Footer**

Via the “Reports” icon you reach your current therapy progress report. Provided you gave your consent during onboarding and complete the questionnaires regularly, the App can create a therapy progress report for you.

On the same page you will also find a downward arrow at the top right: if you tap the arrow, Vertidisan creates a PDF file of your current therapy progress report. Vertidisan also offers you the option of exporting the report. A secure export is possible to your electronic patient record (ePA) and to the practice system of your statutory health insurance physician, provided they allow the import.

You can also export the report via the menu.

Via the “Answer questionnaire” button you reach the questionnaires you selected during onboarding. The questions are put to you in the form of a written dialogue. If possible, complete the questionnaires you have selected once a week.

### **Knowledge Icon in the Footer**

Via the “Knowledge” icon you reach the knowledge library. The knowledge library offers you many specialist medical articles on the subject of vertigo, literature references, quiz questions, the instructions for evaluating your therapy report, Otago exercises and the therapy unit “End of therapy”.

- **Tiles**

By clicking on a tile, you can open the respective content.

- **PDF**

You can also reach the PDF documents via the tiles. Some of the tiles have a downward arrow. If you click on the arrow, you can download the PDF document.

## **Menu Icon in the Footer**

In the menu you will find the following menu items:

- **Account**

- Here you can change your e-mail address and your password.
- Here you can export your data. A secure export is possible to your ePA and to the system of your statutory health insurance physician, if they allow the import. Details can be found in the privacy policy.
- Here you can end your therapy. This also deletes your data. However, you can save your therapy report as a PDF beforehand.
- Here you can delete your data. However, you can save your therapy report as a PDF beforehand.

- **Consents/Withdrawal**

- Notification. Daily reminder function if you have not yet carried out your therapy. For this, you must additionally activate the desired sounds and the desired notification on your device. By removing the tick, you can switch off the reminder.
- Terms of use
- Privacy policy
- Processing of personal data
- Processing of health data. If you have consented, you can **simply withdraw your consent** here by removing the tick. Your health data will then be deleted. Beforehand, you can call them up as a PDF via the "Reports" button in the footer and save them outside the App. Despite the withdrawal, the App will continue to function normally afterwards; only the therapy report will be missing. Details can be found in the privacy policy.

- **Legal**

- Privacy policy. **Read how elaborate and comprehensive Vertidisan's data protection is.**
- General Terms and Conditions
- Legal notice (imprint)
- Information for doctors.
- Side effect reports

- About us.
- **Product information**
  - Product overview. A short introduction to Vertidisan.
  - Areas of application. When and for what Vertidisan is used and what the **contraindications** are.
  - Side effects. Including **side effect reports**
  - Package leaflet. Also referred to as “information for use”.
  - Instructions for use, also user manual
- **Medical device**
  - Name
  - Version. Version number and language (D = German)
  - UDI. For identification as a medical device.
  - Medical device manufacturer
  - Release
  - Intended purpose. You may only use the App for its intended use/intended purpose.
- **Actions**
  - Contact support
  - Log out
- **Quiz questions**

You can reach the quiz questions via the “Knowledge” icon.
- **Evaluation of your therapy report**

You can reach the instructions for evaluating your therapy report via the “Knowledge” icon in the footer. If you hand your therapy report to your doctor, it is advisable to also give him/her a printout of the PDF file of the evaluation instructions. You can reach the therapy report itself via the “Reports” button in the footer.
- **End of therapy**

You can reach the therapy unit “End of therapy” via the “Knowledge” icon.

## 8.6 Intended Use

You may only use the App for its intended use. Vertidisan serves to improve vertigo in peripheral vestibular disorders. Vertidisan also serves to document manually entered health data and to forward this documentation. Vertidisan is not intended to produce or provide information that is used to make decisions for diagnostic or therapeutic purposes.

## **9.0. Classification and Operating System Requirements**

In accordance with Article 51 and Annex VII (Chapter 3) of the Medical Device Regulation (MDR), Vertidisan® is a mobile application in the sense of stand-alone software that has no measuring function and does not make therapeutic decisions. Within the meaning of the MDR (2017), it is therefore a Class I medical device. To use Vertidisan®, a smartphone or tablet with Android operating system (version 12 or higher) or Apple iOS operating system (iOS 18 or higher) is required. A working internet connection is required.

## **10.0. Performance Characteristics**

As the manufacturer, Digitineers GmbH & CoKG guarantees an average online availability of Vertidisan® of 99.5 % per year. Vertidisan® has no functional extensions that can be purchased additionally.

## **11.0 Compatibility Commitments Regarding Supported Platforms**

Vertidisan® runs on an internet-enabled device (smartphone or tablet with Android operating system version 12 or higher), an iPhone or an iPad (iOS operating system version 18 or higher), provided it has internet access sufficient for video streaming.

## **12.0 CE Marking**

Vertidisan® is a CE-marked Class I medical device under the EU MDR (Medical Device Regulation 2017) and at the same time meets the requirements of the German Medical Devices Implementation Act (Medizinproduktegesetz, MPDG).

## **13.0 Data Protection**

Vertidisan® complies with the General Data Protection Regulation (GDPR) and sets high standards for the technical protection of your data. We guarantee that your medical data remain exclusively on secure servers in Germany and that we protect them with the greatest effort. More on data protection, which data are collected and what your rights are can be found in the privacy policy via the menu. The software services integrated into Vertidisan® are listed on the Vertidisan® information homepage ([www.vertidisan.de](http://www.vertidisan.de)).

## 14.0 Product Labelling

In accordance with DIN EN ISO 15223-1:2017-04, the CE marking, manufacturer identification and file identification can be found on the website [www.vertidisan.de](http://www.vertidisan.de) and in the App.

## 15.0 Contact

You can send us an e-mail at any time with questions about content or technical matters! Please contact us at [support@digitineers.de](mailto:support@digitineers.de)

## 16.0 Manufacturer

Vertidisan® is manufactured and placed on the market by Digitineers GmbH & CoKG in Germany.

Digitineers GmbH & CoKG  
August-Bebel-Str. 9  
72072 Tübingen  
Germany  
[www.digitineers.com](http://www.digitineers.com)

## Annex 1:

### Indications

#### Vertigo Due to Disorders of Vestibular Function

ICD-10

H81.0 Ménière's disease

H81.1 Benign paroxysmal vertigo

H81.2 Vestibular neuronitis

H81.3 Other peripheral vertigo

Incl.:

Vertigo: aural

Vertigo: otogenic

Vertigo: peripheral NOS

H81.8 Other disorders of vestibular function

H81.9 Disorder of vestibular function, unspecified

Incl.:

Vertiginous syndrome NOS

R29.6 Tendency to fall, not elsewhere classified

Incl.: Tendency to fall due to other unclear health conditions

Tendency to fall in older persons


### Contraindications

#### **Vertidisan® must not be used in the case of**

- Diseases of the cervical spine, cervical spine syndrome
- Acute inflammatory diseases
- Febrile condition
- Heart failure (cardiovascular decompensation)
- Lung failure (pulmonary insufficiency)
- Motor disorders
- Seizure disorders (e.g. epilepsy)
- Spinal cord damage
- Psychoses
- Skin sensitivity disorders with walking impairment
- Thrombosis

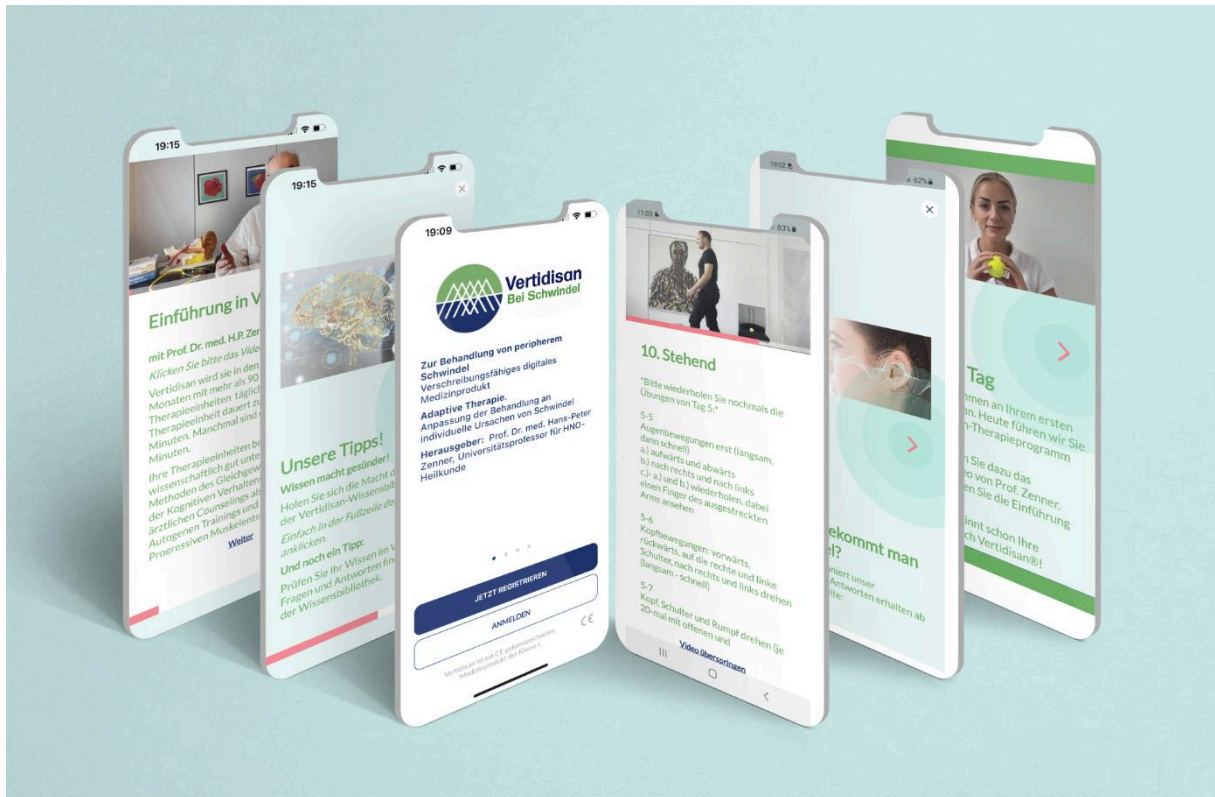- Blood coagulation disorders, including anticoagulant therapy
- Large skin defects
- Burns requiring hospitalisation
- Metastases of malignant tumours
- Severe visual impairment
- Frailty, frailty syndrome
- Pregnancy

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© Copyright 2023 - Digitineers GmbH & Co KG, Tübingen